

ALCOTEST CHECKLIST

Municipality: FLORENCE TWP. Alcotest Ser. #: ARTN - 0053

County: BURLINGTON Date of Calibration: 01-05-2011

1. **Certificate of Accuracy Alcotest 7110 MKIII-C** from Draeger Safety for instrument used in the A.I.R.
2. **Certificate of Accuracy CU34 Unit** on Alcotest Instrument used.
Ser. #: DDUM 53-0222
3. **Certificate of Accuracy Alcotest 7110 Temperature Probe** from Draeger Safety for instrument in A.I.R. or equivalent.
Ser. #: DDWT P2-293

4. **Digital Temperature Measuring System Report of Calibration.**
Ser. #: 101 733 534

5. A. Alcotest 7110 Calibration Record
 B. Alcotest 7110 Calibration Certificate **Part I - Control Tests.**
 C. Alcotest 7110 Calibration Certificate **Part II - Linearity Tests.**
 D. **Alcotest Card** of operator/coordinator who performed tests.
 E. **Certificate of Accuracy Alcotest 7110 Temperature Probe** from Draeger Safety used in the Calibration Tests ["Black Key" probe of Breath Test Coordinator].
Ser. #: DDWA P2-016

6. **Certificates of Analysis** for each **Simulator Solution** used in Calibration/Linearity Tests:
- | | | | | |
|--|-----------------|-----------|----------|------------|
| <input checked="" type="checkbox"/> A. | 0.04% Solution. | <u>10</u> | <u>A</u> | <u>073</u> |
| <input checked="" type="checkbox"/> B. | 0.08% Solution. | <u>10</u> | <u>A</u> | <u>074</u> |
| <input checked="" type="checkbox"/> C. | 0.10% Solution. | <u>10</u> | <u>F</u> | <u>080</u> |
| <input checked="" type="checkbox"/> D. | 0.16% Solution. | <u>10</u> | <u>A</u> | <u>075</u> |

7. Certificate of Accuracy Alcotest CU34 Simulators from Draeger Safety (when conducting the Calibration/Linearity Tests) for:
- A. 0.04% used in Calibration/Linearity Testing. DDWE 53-0196
- B. 0.08% used in Calibration/Linearity Testing. DDWE 53-0203
- C. 0.10% used in Calibration/Linearity Testing. DDUM 53-0222
[Same as instrument's CU34 unit].
- D. 0.16% used in Calibration/Linearity Testing. DDWE 53-0205

8. A. New Standard Solution Report following Calibration.
 B. Calibrating CU34 Unit for same [same as CU34 unit on instrument].
 C. Certificate of Analysis 0.10% solution for same.
Lot #: 10 F 079
 D. Alcotest card of operator/coordinator who completed change.

Dräger safety

Alcotest[®] 7110 MKIII-C

CERTIFICATE OF ACCURACY

This is to certify that the Alcotest 7110 MKIII-C has been tested for accuracy and found to be in compliance with the National Highway Traffic Safety Administration Standard for evidential breath testing devices. The Alcotest MKIII-C is compliant as a "mobile" and "nonmobile" EBT with 49 FR 48884, 49 FR 48884 and 58 FR 48705. The manufacturer recommends accuracy verification of this instrument within 12 months of the calibration date below, or sooner, according to your State Specifications.

Certification Date:

11-18-2004

SERIAL NUMBER:

ARTN-0053

Draeger Safety Diagnostics, Inc.
Durango, CO





Drägersafety

CERTIFICATE OF ACCURACY

This Certificate of Accuracy verifies that the specified unit has been examined and found to be in compliance with National Highway and Traffic Safety Administration regulations for devices used to calibrate Evidential Breath Testers.
(F.R. Vol. 59 No. 249 12/19/94 Notices)

Draeger Safety Diagnostics, Inc.

- Model: ALCOTEST® CU34
- Model: MARK IIA
- Other: _____

Serial Number:

DDUMS3-0222

Certification Date

11-17-10

Technician

mm

Re-Certification Due Date

11-17-11

Drägersafety

ALCOTEST® 7110 TEMPERATURE PROBE

CERTIFICATE OF ACCURACY

This is to certify that the Alcotest® 7110 Temperature Probe has been tested for accuracy with instrumentation that is traceable to the National Institute of Standards and Technology (NIST).

The manufacturer recommends accuracy verification of the Temperature Probe within 12 months of the certification date below, or sooner, according to your State Specification.

For accurate temperature readings, the probe value on this certificate, noted below, must be programmed into the Alcotest® 7110.

Serial Number Temp. Probe

DDWTP2-293

Certification date:

11-17-10

Next Certification due:

11-17-11

Probe Value

105

Draeger Safety Diagnostics, Inc.
Technical Service Department

mm



Calibration
Certificate No. 1750.01

Calibration complies with ISO 9001
ISO/IEC 17025 AND ANSI/NCSL Z540-1



Cert. No.: 4000-2966289

Traceable® Certificate of Calibration for Digital Thermometer

Manufactured for and distributed by: VWR International, P.O. Box 2158, Secaucus, NJ 07094 U.S.A.

Instrument Identification:

Model: 61220-601

S/N: 101733534

Manufacturer : Control Company

Standards/Equipment:

Description	Serial Number	Due Date	NIST Traceable Reference
Temperature Calibration Bath TC-179	A45240		
Thermistor Module	A17118	11/19/10	A9B21010
Temperature Probe	128	12/10/10	A9B23079
Temperature Calibration Bath TC-231	A79341		
Temperature Probe	3039	12/10/10	A9B23080-1
Temperature Calibration Bath TC-218	A73332		
Thermistor Module	A27129	7/09/10	1000264338
Temperature Probe	5202	3/11/11	B0310050
Temperature Calibration Bath TC-256	B01375		
Temperature Probe	157	7/27/10	A9708011-4

Certificate Information:

Technician: 104

Procedure: CAL-06

Cal Date: 6/08/10

Cal Due: 6/08/12

Test Conditions: 24.0°C 42.0 %RH 1015 mBar

Calibration Data: (New Instrument)

Unit(s)	Nominal	As Found	In Tol	Nominal	As Left	In Tol	Min	Max	±U	TUR
°C		N.A.		0.002	-0.001	Y	-0.048	0.052	0.013	3.8:1
°C		N.A.		24.999	25.000	Y	24.949	25.049	0.013	3.8:1
°C		N.A.		60.001	59.998	Y	59.951	60.051	0.018	2.8:1
°C		N.A.		100.001	100.002	Y	99.951	100.051	0.013	3.8:1

This instrument was calibrated using instruments traceable to National Institute of Standards and Technology.

A Test Uncertainty Ratio of at least 4:1 is maintained unless otherwise stated and is calculated using the expanded measurement uncertainty. Uncertainty evaluation includes the instrument under test and is calculated in accordance with the ISO "Guide to the Expression of Uncertainty in Measurement" (GUM). The uncertainty represents an expanded uncertainty using a coverage factor k=2 to approximate a 95% confidence level. In tolerance conditions are based on test results falling within specified limits with no reduction by the uncertainty of the measurement. The results contained herein relate only to the item calibrated. This certificate shall not be reproduced except in full, without written approval of Control Company.

Nominal=Standard's Reading; As Left=Instrument's Reading; In Tol=In Tolerance; Min/Max=Acceptance Range; ±U=Expanded Measurement Uncertainty; TUR=Test Uncertainty Ratio; Accuracy=(Max-Min)/2; Min = Nominal(Rounded) - Tolerance; Max = Nominal(Rounded) + Tolerance. Date=MM/DD/YY

Nicol Rodriguez
Nicol Rodriguez, Quality Manager

Wallace Berry
Wallace Berry, Technical Manager

Maintaining Accuracy:

In our opinion once calibrated your Digital Thermometer should maintain its accuracy. There is no exact way to determine how long calibration will be maintained. Digital Thermometers change little, if any at all, but can be affected by aging, temperature, shock, and contamination.

Recalibration:

For factory calibration and re-certification traceable to National Institute of Standards and Technology contact Control Company.

CONTROL COMPANY 4455 Rex Road Friendswood, TX 77546 USA
Phone 281 482-1714 Fax 281 482-9448 service@control3.com www.control3.com

Control Company is an ISO 17025:2005 Calibration Laboratory Accredited by (A2LA) American Association for Laboratory Accreditation, Certificate No. 1750.01.
Control Company is ISO 9001:2008 Quality Certified by (DNV) Det Norske Veritas, Certificate No. CERT-01805-2008-AQ-HOU-ANAB.
International Laboratory Accreditation Cooperation (ILAC) - Multilateral Recognition Arrangement (MRA).

Alcotest 7110 Calibration Record

Equipment Alcotest 7110 MKIII-C Serial No.: AR1N-0053
Location: FLORENCE TOWNSHIP P.D.
Calibration File No.: 00427 **Calib. Date:** 01/05/2011 **Calib. No.:** 00010
Certification File No.: 00393 **Cert. Date:** 07/06/2010 **Cert. No.:** 00008
Linearity File No.: 00394 **Lin. Date:** 07/06/2010 **Lin. No.:** 00008
Solution File No.: 00424 **Soln. Date:** 12/22/2010 **Soln. No.:** 00096
Sequential File No.: 00427 **File Date:** 01/05/2011

Calibrating Unit: WET **Model No.:** CU-34 **Serial No.:** DDUM S3-0222
Control Solution %: 0.100% **Expires:** 06/14/2012
Solution Control Lot: 10F080 **Bottle No.:** 0381

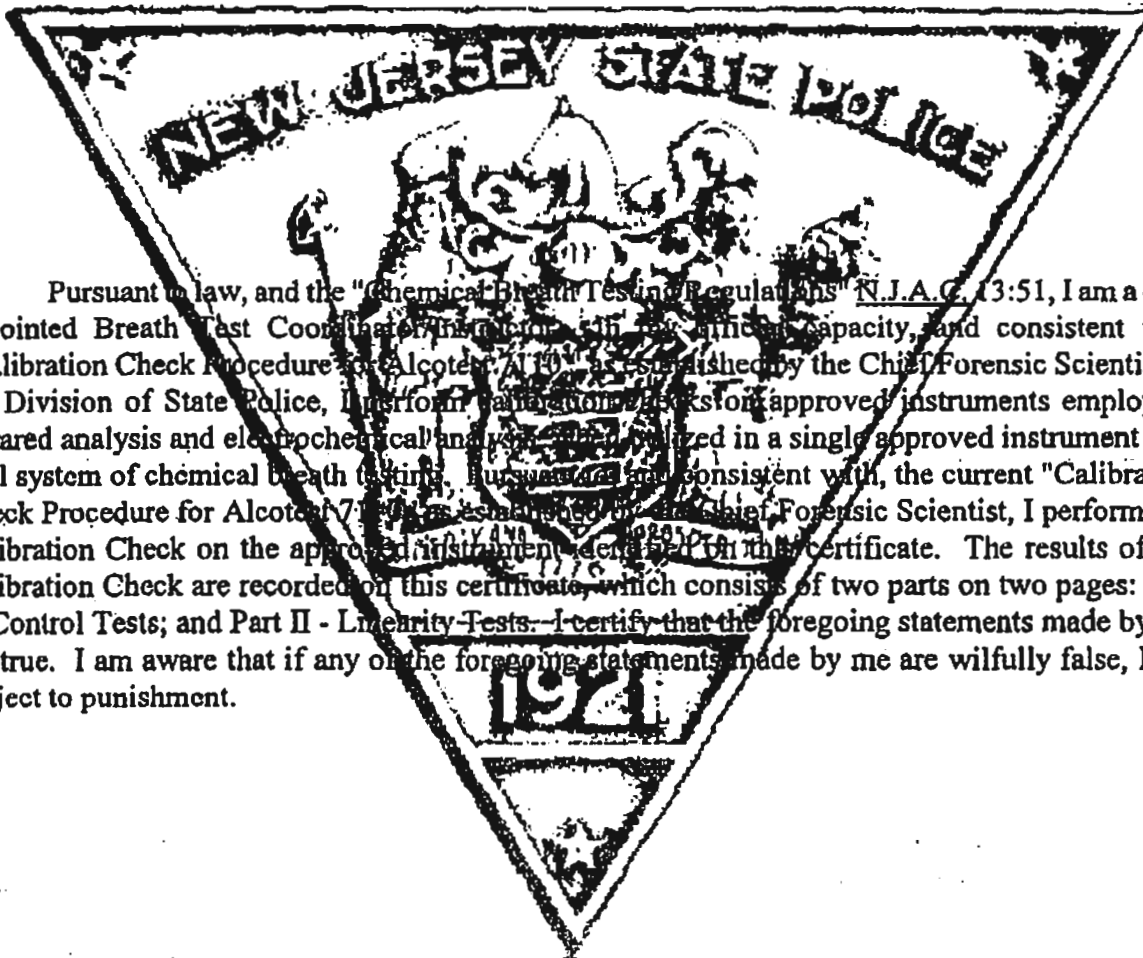
Coordinator

Last Name: DELLANOCE **First Name:** JOSEPH **MI:** S

Signature: *J. Dellanoce # 6027* **Badge No.:** 6027
Date: 01/05/2011

*Black Key Temperature Probe Serial.....# DDWAP2-016 *jo*

*Digital NIST Temperature Measuring System Serial.....# 101733534 *jo*



Pursuant to law, and the "Chemical Breath Testing Regulations" N.J.A.C. 13:51, I am a duly appointed Breath Test Coordinator in my official capacity, and consistent with "Calibration Check Procedure for Alcotest 7110" as established by the Chief Forensic Scientist of the Division of State Police, I perform calibration checks on approved instruments employing infrared analysis and electrochemical analysis, which are analyzed in a single approved instrument as a dual system of chemical breath testing. Pursuant and consistent with, the current "Calibration Check Procedure for Alcotest 7110" as established by the Chief Forensic Scientist, I performed a Calibration Check on the approved instrument identified on this certificate. The results of my Calibration Check are recorded on this certificate, which consists of two parts on two pages: Part I - Control Tests; and Part II - Linearity Tests. I certify that the foregoing statements made by me are true. I am aware that if any of the foregoing statements made by me are wilfully false, I am subject to punishment.

Alcotest 7110 Calibration Certificate

Part I - Control Tests

Equipment Alcotest 7110 MKIII-C Serial No.: ARTN-0053
Location: FLORENCE TOWNSHIP P.D.
Calibration File No.: 00427 **Calib. Date:** 01/05/2011 **Calib. No.:** 00010
Certification File No.: 00428 **Cert. Date:** 01/05/2011 **Cert. No.:** 00009
Linearity File No.: 00394 **Lin. Date:** 07/06/2010 **Lin. No.:** 00008
Solution File No.: 00424 **Soln. Date:** 12/22/2010 **Soln. No.:** 00096
Sequential File No.: 00428 **File Date:** 01/05/2011

Calibrating Unit: WET **Model No.:** CU-34 **Serial No.:** DDUM S3-0222
Control Solution %: 0.100% **Expires:** 06/14/2012
Solution Control Lot: 10F080 **Bottle No.:** 0381

Function	Result %BAC	Time HH:MM	Temperature Simulator (°C)	Comment(s) or Error(s)
Ambient Air Blank	0.000%	11:02S		
Control 1 EC	0.099%	11:02S	34.0°C	*** TEST PASSED ***
Control 1 IR	0.100%	11:02S	34.0°C	*** TEST PASSED ***
Ambient Air Blank	0.000%	11:03S		
Control 2 EC	0.098%	11:04S	34.0°C	*** TEST PASSED ***
Control 2 IR	0.099%	11:04S	34.0°C	*** TEST PASSED ***
Ambient Air Blank	0.000%	11:05S		
Control 3 EC	0.098%	11:05S	34.0°C	*** TEST PASSED ***
Control 3 IR	0.098%	11:05S	34.0°C	*** TEST PASSED ***
Ambient Air Blank	0.000%	11:06S		

All tests within acceptable tolerance.

Coordinator
Last Name: DELLANOCE First Name: JOSEPH
Signature: *[Handwritten Signature]* Badge No. 1115
Date: 01/05/2011

Pursuant to law, and the "Chemical Breath Testing Regulations" N.J.A.C. 13:51, I am a duly appointed Breath Test Coordinator in the above capacity, and consistent with "Calibration Check Procedure for Alcotest 7110" as established by the Chief Forensic Scientist of the Division of State Police, I performed calibration checks on approved instruments employing infrared analysis and electrochemical analysis which are utilized in a single approved instrument as a dual system of chemical breath testing. Pursuant and consistent with, the current "Calibration Check Procedure for Alcotest 7110" as established by the Chief Forensic Scientist, I performed a Calibration Check on the approved instrument identified on this certificate. The results of my Calibration Check are recorded on this certificate, which consists of two parts on two pages: Part I - Control Tests; and Part II - Linearity Tests. I certify that the foregoing statements made by me are true. I am aware that if any of the foregoing statements made by me are wilfully false, I am subject to punishment.

Alcotest 7110 Calibration Certificate

Part II - Linearity Tests

Equipment Alcotest 7110 MKIII-C Serial No.: ARTN-0053
Location: FLORENCE TOWNSHIP P.D.
Calibration File No.: 00427 **Calib. Date:** 01/05/2011 **Calib. No.:** 00010
Certification File No.: 00428 **Cert. Date:** 01/05/2011 **Cert. No.:** 00009
Linearity File No.: 00429 **Lin. Date:** 01/05/2011 **Lin. No.:** 00009
Solution File No.: 00424 **Soln. Date:** 12/22/2010 **Soln. No.:** 00096
Sequential File No.: 00429 **File Date:** 01/05/2011

Calibrating Unit: WET **Model No.:** CU-34 **Serial No.:** DDWE S3-0196
Control Solution %: 0.040% **Expires:** 01/12/2012
Solution Control Lot: 10A073 **Bottle No.:** 1127

Calibrating Unit: WET **Model No.:** CU-34 **Serial No.:** DDWE S3-0203
Control Solution %: 0.080% **Expires:** 01/15/2012
Solution Control Lot: 10A074 **Bottle No.:** 1011

Calibrating Unit: WET **Model No.:** CU-34 **Serial No.:** DDWE S3-0205
Control Solution %: 0.160% **Expires:** 01/21/2012
Solution Control Lot: 10A075 **Bottle No.:** 0591

Function	Result %BAC	Time HH:MM	Temperature Simulator (°C)	Comment(s) or Error(s)
Ambient Air Blank	0.000%	11:23S		
Control 1 EC	0.041%	11:23S	34.0°C	*** TEST PASSED ***
Control 1 IR	0.040%	11:23S	34.0°C	*** TEST PASSED ***
Ambient Air Blank	0.000%	11:25S		
Control 2 EC	0.042%	11:26S	34.0°C	*** TEST PASSED ***
Control 2 IR	0.040%	11:26S	34.0°C	*** TEST PASSED ***
Ambient Air Blank	0.000%	11:27S		
Control 3 EC	0.083%	11:28S	34.0°C	*** TEST PASSED ***
Control 3 IR	0.078%	11:28S	34.0°C	*** TEST PASSED ***
Ambient Air Blank	0.000%	11:30S		
Control 4 EC	0.080%	11:30S	34.0°C	*** TEST PASSED ***
Control 4 IR	0.078%	11:30S	34.0°C	*** TEST PASSED ***
Ambient Air Blank	0.000%	11:32S		
Control 5 EC	0.161%	11:32S	34.0°C	*** TEST PASSED ***
Control 5 IR	0.156%	11:32S	34.0°C	*** TEST PASSED ***
Ambient Air Blank	0.000%	11:34S		
Control 6 EC	0.158%	11:35S	34.0°C	*** TEST PASSED ***
Control 6 IR	0.157%	11:35S	34.0°C	*** TEST PASSED ***
Ambient Air Blank	0.000%	11:37S		

All tests within acceptable tolerance.

Coordinator

Last Name: DELLANOCE

First Name: JOSEPH

M1: S

Signature: _____

Tr II J. Dellanoce # 6027

Badge No.: 6027

Date: 01/05/2011

DEPARTMENT OF
Justice and Public Safety
 This is to certify that

Joseph S. Dellanoce
 New Jersey State Police

IS QUALIFIED AND COMPETENT TO CONDUCT CERTIFIED ALCOHOL ANALYSES PURSUANT TO CHAPTER 143 OF
 TITLE 14 OF THE STATUTES OF THE STATE OF NEW JERSEY IN ACCORDANCE WITH THE PROVISIONS OF
 A METHOD TO DETERMINE INTOXICATION
 OF PERSONS UNDER THE INFLUENCE OF ALCOHOL BY THE METHOD OF BLOW INTO TEST 7110 MKIII-C

TWO THOUSAND AND FIVE

25th DAY OF March

[Signature]
 SUPERVISOR
 NEW JERSEY STATE POLICE

[Signature]
 ATTORNEY GENERAL
 STATE OF NEW JERSEY

ORIGINAL COURSE DATES

DATE	Refresher Course PLACE	INSTRUCTOR
1. 2/15/07	CCPA	C. Kozell
2. 12/7/09	DLPA	William S. ...
3.		
4.		
5.		
6.		
7.		
8.		
9.		

S.P. 2378 (Rev. 12/04)

DEPARTMENT OF
Fire and Public Safety
 While it is certified that
Joseph D. Dallanocce,
 New Jersey State Police

IS QUALIFIED AND CAPABLE TO PERFORM THE ANALYSIS OF EVIDENCE IN THE
 FIELD OF FORENSIC SCIENCE BY USING THE **Analyzer**
 AND TO DETECT THE PRESENCE OF

ON THE 13th DAY OF Dec.
 TWO THOUSAND AND 01

C. P. Felt
 SUPERVISOR
 NEW JERSEY STATE POLICE

J. P. J.
 ATTORNEY GENERAL
 STATE OF NEW JERSEY

ORIGINAL COURSE DATES: 8/90 - 9/3/99

Refresher Course

DATE	PLACE	INSTRUCTOR
1. 12/12/01	A. G. C.	ROTTBE.
2. 2-20-03	Field	P. J. J.
3. 10-3-03	Field	C. J. J.
4.		
5.		
6.		
7.		
8.		
9.		

SP-297B (Rev. 10/99)

DEPARTMENT OF
Traffic and Public Safety
It is to certify that

Joseph S. Dellanoce
New Jersey State Police

IS QUALIFIED AND COMPETENT TO CONDUCT CHEMICAL ANALYSIS AND THIS PURSUANT TO CHAPTER 80 OF
THE LAWS OF 1966 IN THE CREATION OF THE Breath Test Coordinator/Instructor
A METHOD TO DETERMINE INTOXICATION
GIVEN UNDER MY HAND AT TRENTON, NEW JERSEY, THIS 14th DAY OF January

TWO THOUSAND AND EIGHT


SUPERINTENDENT
NEW JERSEY STATE POLICE


ATTORNEY GENERAL
STATE OF NEW JERSEY

ORIGINAL COURSE DATES

DATE	Refresher Course PLACE	INSTRUCTOR
1. <u>12/7/09</u>	<u>OCPA</u>	<u>Adam Standa</u>
2. _____	_____	_____
3. _____	_____	_____
4. _____	_____	_____
5. _____	_____	_____
6. _____	_____	_____
7. _____	_____	_____
8. _____	_____	_____
9. _____	_____	_____



State of New Jersey
OFFICE OF THE ATTORNEY GENERAL
DEPARTMENT OF LAW AND PUBLIC SAFETY
PO BOX 080
TRENTON, NJ 08625-0080

JON S. CORZINE
Governor

ANNE MILGRAM
Attorney General

January 14, 2008

Col. Joseph R. Fuentes, Superintendent
Division of State Police
Division Headquarters
P.O. Box 7068
West Trenton, New Jersey 08628

Re: Breath Test Coordinator/Instructor, Certification -
Trooper Joseph S. Dellanoce #6027

Dear Colonel Fuentes:

Pursuant to the provisions of *N.J.A.C. 13:51-2.1 (b) and (c)*, as adopted and promulgated under the provisions of *N.J.S.A. 39:4-50.3, 39:3-10.25 and 12:7-56*, I hereby approve Trooper Joseph S. Dellanoce #6027, as a duly certified Breath Test Coordinator/Instructor. This approval is effective immediately.

Very truly yours,

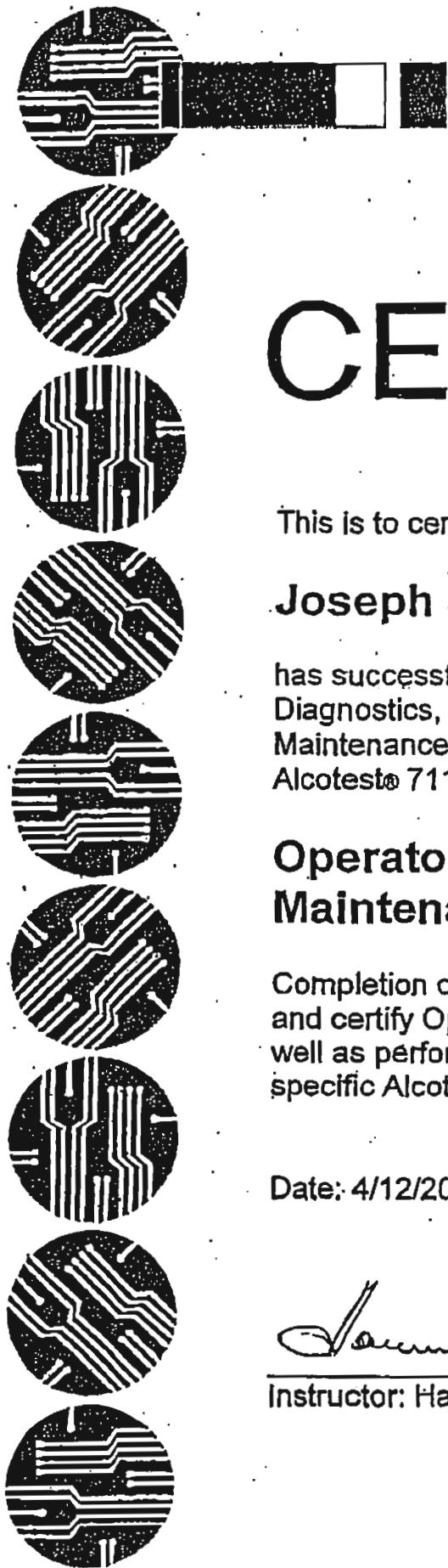

Anne Milgram
Attorney General

jk

c. Trooper Joseph S. Dellanoce #6027, Alcohol/Drug Test Unit, Division of State Police
Lt. Mark Kolodzieski, Unit Head, Alcohol/Drug Test Unit, Division of State Police



329474



CERTIFICATE

This is to certify that

Joseph S. Dellanoce #6027

has successfully completed the two day Draeger Safety Diagnostics, Inc. Operator Training and Preventive Maintenance Course on the New Jersey specific Alcotest® 7110 MKIII-C and is hereby certified as a qualified

**Operator Trainer and
Maintenance Technician**

Completion of this course qualifies this individual to train and certify Operators in the proper use and operation as well as perform Preventive Maintenance on the New Jersey specific Alcotest® 7110 MKIII-C.

Date: 4/12/2007

Instructor: Hansuelli Rysler

Dräger

Dräger safety

ALCOTEST® 7110 TEMPERATURE PROBE

CERTIFICATE OF ACCURACY

This is to certify that the Alcotest® 7110 Temperature Probe has been tested for accuracy with instrumentation that is traceable to the National Institute of Standards and Technology (NIST). The manufacturer recommends accuracy verification of the Temperature Probe within 12 months of the certification date below, or sooner, according to your State Specification. For accurate temperature readings, the probe value on this certificate, noted below, must be programmed into the Alcotest® 7110.

Serial Number Temp. Probe

DDWAP2-016

Certification date:

2-3-10

Next Certification due:

2-3-11

Probe Value

106

Dräger Safety Diagnostics, Inc.
Technical Service Department

mm



State of New Jersey
 OFFICE OF THE ATTORNEY GENERAL
 DEPARTMENT OF LAW AND PUBLIC SAFETY
 DIVISION OF STATE POLICE
 Post Office Box 7068
 WEST TRENTON NJ 08628-0068
 (609) 882-2000

PAULA T. DOW
 Acting Attorney General

COLONEL JOSEPH R. FUENTES
 Superintendent

CHRIS CHRISTIE
 Governor

CERTIFICATION OF ANALYSIS
0.04 PERCENT BREATH ALCOHOL SIMULATOR SOLUTION

ACCEPTANCE SPECIFICATIONS FOR BREATH ALCOHOL SIMULATOR SOLUTION: Ethyl alcohol concentration within, but not exceeding, the range of 0.045 to 0.051 grams per 100 milliliters of solution.


MANUFACTURER: Drager Safety, Inc. ANALYSIS DATE: 2/2/2010

BREATH ALCOHOL SIMULATOR SOLUTION LOT NUMBER: 10A073

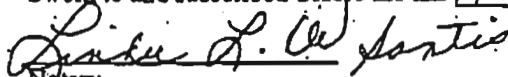
Representative samples of the above-referenced Lot Number were tested by Gas Chromatography and found to have a mean ethyl alcohol concentration range of 0.0479 to 0.0481 grams per 100 milliliters of solution.

This lot of breath alcohol simulator solution may be utilized as a known traceable standard for the purpose of conducting periodic tests, pursuant to N.J.A.C. 13:51-3.4, of approved breath test instruments (N.J.A.C. 13:51-3.5) utilized by law enforcement agencies in this State. The manufacturer's expiration date for this lot of breath alcohol simulator solution is January 12, 2012.

As Assistant Chief Forensic Scientist of the Division of State Police, I hereby certify and attest that the tests and results documented in this Certificate of Analysis were performed at my direction and under my supervision by personnel of, and at, the Office of Forensic Sciences of the Division of State Police on properly functioning and calibrated instruments and equipment. All procedures utilized are accurate, objective, and performed on a routine basis by personnel within the Office of Forensic Sciences, in accordance with their professional duties and responsibilities.


 Kenneth W. Krawalek, M.S.
 Assistant Chief Forensic Scientist
 Division of State Police

Sworn to and subscribed before me this 19th day of February, 2010.


 Notary

Linda L. DeBartis
 Notary Public, New Jersey
 My Commission Expires 8-17-14



New Jersey Is An Equal Opportunity Employer
 Printed on Recycled Paper and Recyclable





State of New Jersey
 OFFICE OF THE ATTORNEY GENERAL
 DEPARTMENT OF LAW AND PUBLIC SAFETY
 DIVISION OF STATE POLICE
 Post Office Box 7068
 West Trenton NJ 08628-0068
 (609) 882-2000

CHRIS CHRISTIE
 Governor

PAULA T. DOW
 Acting Attorney General

COLONEL JOSEPH R. FUENTES
 Superintendent

CERTIFICATION OF ANALYSIS
0.08 PERCENT BREATH ALCOHOL SIMULATOR SOLUTION

ACCEPTANCE SPECIFICATIONS FOR BREATH ALCOHOL SIMULATOR SOLUTION: Ethyl alcohol concentration within, but not exceeding, the range of 0.094 to 0.099 grams per 100 milliliters of solution.

MANUFACTURER: Drager Safety, Inc.

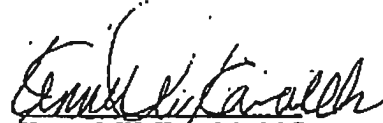
ANALYSIS DATE: 2/3/2010

BREATH ALCOHOL SIMULATOR SOLUTION LOT NUMBER: 10A074

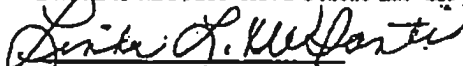
Representative samples of the above-referenced Lot Number were tested by Gas Chromatography and found to have a mean ethyl alcohol concentration range of 0.0954 to 0.0958, grams per 100 milliliters of solution.

This lot of breath alcohol simulator solution may be utilized as a known traceable standard for the purpose of conducting periodic tests, pursuant to N.J.A.C. 13:51-3.4, of approved breath test instruments (N.J.A.C. 13:51-3.5) utilized by law enforcement agencies in this State. The manufacturer's expiration date for this lot of breath alcohol simulator solution is January 15, 2012.

As Assistant Chief Forensic Scientist of the Division of State Police, I hereby certify and attest that the tests and results documented in this Certificate of Analysis were performed at my direction and under my supervision by personnel of, and at, the Office of Forensic Sciences of the Division of State Police on properly functioning and calibrated instruments and equipment. All procedures utilized are accurate, objective, and performed on a routine basis by personnel within the Office of Forensic Sciences, in accordance with their professional duties and responsibilities.


 Kenneth W. Kawalek, M.S.
 Assistant Chief Forensic Scientist
 Division of State Police

Sworn to and subscribed before me this 19th day of February, 2010.


 Notary.

Linda L. DeSantis
 Notary Public, New Jersey
 My Commission Expires 8-17-14



New Jersey Is An Equal Opportunity Employer
 Printed on Recycled Paper and Recyclable





State of New Jersey
 OFFICE OF THE ATTORNEY GENERAL
 DEPARTMENT OF LAW AND PUBLIC SAFETY
 DIVISION OF STATE POLICE
 POST OFFICE BOX 7068
 WEST TRIGNION NJ 08628-0068
 (609) 882-2000

CHRIS CHRISTIE
 Governor

KIM GUADAGNO
 Lt. Governor

PAULA T. DOW
 Attorney General

COLONEL JOSEPH R. FUENTES
 Superintendent

CERTIFICATION OF ANALYSIS
0.10 PERCENT BREATH ALCOHOL SIMULATOR SOLUTION

ACCEPTANCE SPECIFICATIONS FOR BREATH ALCOHOL SIMULATOR SOLUTION: Ethyl alcohol concentration within, but not exceeding, the range of 0.117 to 0.125 grams per 100 milliliters of solution.

MANUFACTURER: Draeger Safety, Inc.

ANALYSIS DATE: 7/8/2010

BREATH ALCOHOL SIMULATOR SOLUTION LOT NUMBER: 10F080

Representative samples of the above-referenced Lot Number were tested by Gas Chromatography and found to have a mean ethyl alcohol concentration range of 0.1197 to 0.1267 grams per 100 milliliters of solution.

This lot of breath alcohol simulator solution may be utilized as a known traceable standard for the purpose of conducting periodic tests, pursuant to N.J.A.C. 13:51-3.4, of approved breath test instruments (N.J.A.C. 13:51-3.5) utilized by law enforcement agencies in this State. The manufacturer's expiration date for this lot of breath alcohol simulator solution is June 14, 2012.

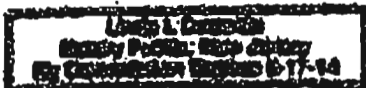
As Chief Forensic Scientist of the Division of State Police, I hereby certify and attest that the tests and results documented in this Certificate of Analysis were performed at my direction and under my supervision by personnel of, and at, the Office of Forensic Sciences of the Division of State Police on properly functioning and calibrated instruments and equipment. All procedures utilized are accurate, objective, and performed on a routine basis by personnel within the Office of Forensic Sciences, in accordance with their professional duties and responsibilities.

Ajit R. Tungare

Ajit R. Tungare
 Chief Forensic Scientist
 Division of State Police

Sworn to and subscribed before me this 5th day of July, 2010.

Linda L. DeFuentes
 Notary



New Jersey Is An Equal Opportunity Employer
 Printed on Recycled Paper and Recyclable





CHRIS CHRISTIE
Governor

State of New Jersey
OFFICE OF THE ATTORNEY GENERAL
DEPARTMENT OF LAW AND PUBLIC SAFETY
DIVISION OF STATE POLICE
POST OFFICE BOX 7068
WEST TRENTON NJ 08628-0068
(609) 882-2000

PAULA T. DOW
Acting Attorney General

COLONEL JOSEPH R. FUENTES
Superintendent

CERTIFICATION OF ANALYSIS
0.16 PERCENT BREATH ALCOHOL SIMULATOR SOLUTION

ACCEPTANCE SPECIFICATIONS FOR BREATH ALCOHOL SIMULATOR SOLUTION: Ethyl alcohol concentration within, but not exceeding, the range of 0.188 to 0.199 grams per 100 milliliters of solution.

MANUFACTURER: Dräger Safety, Inc.


ANALYSIS DATE: 2/4/2010

BREATH ALCOHOL SIMULATOR SOLUTION LOT NUMBER: 10A075

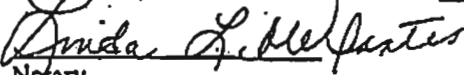
Representative samples of the above-referenced Lot Number were tested by Gas Chromatography and found to have a mean ethyl alcohol concentration range of 0.1913 to 0.1919 grams per 100 milliliters of solution.

This lot of breath alcohol simulator solution may be utilized as a known traceable standard for the purpose of conducting periodic tests, pursuant to N.J.A.C. 13:51-3.4, of approved breath test instruments (N.J.A.C. 13:51-3.5) utilized by law enforcement agencies in this State. The manufacturer's expiration date for this lot of breath alcohol simulator solution is January 21, 2012.

As Assistant Chief Forensic Scientist of the Division of State Police, I hereby certify and attest that the tests and results documented in this Certificate of Analysis were performed at my direction and under my supervision by personnel of, and at, the Office of Forensic Sciences of the Division of State Police on properly functioning and calibrated instruments and equipment. All procedures utilized are accurate, objective, and performed on a routine basis by personnel within the Office of Forensic Sciences, in accordance with their professional duties and responsibilities.


Kenneth W. Kawalek, M.S.
Assistant Chief Forensic Scientist
Division of State Police

Sworn to and subscribed before me this 19th day of February, 2010.


Notary

Linda L. Dondria
Notary Public, New Jersey
My Commission Expires 8-17-14



New Jersey Is An Equal Opportunity Employer
Detailed on Separated Form and Beneficial





Dräger safety

CERTIFICATE OF ACCURACY

This Certificate of Accuracy verifies that the specified unit has been examined and found to be in compliance with National Highway and Traffic Safety Administration regulations for devices used to calibrate Evidential Breath Testers.

(P.R. Vol. 59 No. 249-12/19/94 Notices)

Draeger Safety Diagnostics, Inc.

Model: ALCOTEST® CU34

Model: MARK IIA

Other: _____

Serial Number:

DD W E S 3 - 0196

Certification Date

Technician

Re-Certification Due Date

2/3/10

DM

2/3/11



Draeger safety

CERTIFICATE OF ACCURACY

This Certificate of Accuracy verifies that the specified unit has been examined and found to be in compliance with National Highway and Traffic Safety Administration regulations for devices used to calibrate Evidential Breath Testers.
(F.R. Vol. 59 No. 249 12/19/94 Notices)

Draeger Safety Diagnostics, Inc.

- Model: ALCOTEST® CU34
- Model: MARK II A
- Other: _____

Serial Number:

DDWES3-0203

Certification Date

Technician

Re-Certification Due Date

2/3/10

DM

2/3/11



Draeger safety

CERTIFICATE OF ACCURACY

This Certificate of Accuracy verifies that the specified unit has been examined and found to be in compliance with National Highway and Traffic Safety Administration regulations for devices used to calibrate Evidential Breath Testers.

(F.R. Vol. 59 No. 249 12/19/94 Notices)

Draeger Safety Diagnostics, Inc.

Model: ALCOTEST® CU34

Model: MARK IIA

Other: _____

Serial Number:

DDWES3-0205

Certification Date

Technician

Re-Certification Due Date

2/3/10

DM

2/3/11

Calibrating Unit New Standard Solution Report

Equipment	Alcotest 7110 MKIII-C	Serial No.:	ARTN-0053
Location:	FLORENCE TOWNSHIP P.D.	Calib. No.:	00010
Calibration File No.:	00427	Calib. Date:	01/05/2011
Certification File No.:	00428	Cert. No.:	00009
Linearity File No.:	00429	Lin. Date:	01/05/2011
Solution File No.:	00430	Soln. No.:	00097
Sequential File No.:	00430	File Date:	01/05/2011

Calibrating Unit:	WET	Model No.:	CU-34	Serial No.:	DDUM S3-0222
Control Solution %:	0.100%	Expires:	06/10/2012	Bottle No.:	1094
Solution Control Lot:	10F079				

Function	Result %BAC	Time HH:MM	Temperature Simulator (°C)	Comment(s) or Error(s)
Ambient Air Blank	0.000%	12:44S		
Control 1 EC	0.100%	12:44S	34.0°C	*** TEST PASSED ***
Control 1 IR	0.099%	12:44S	34.0°C	*** TEST PASSED ***
Ambient Air Blank	0.000%	12:45S		
Control 2 EC	0.098%	12:46S	34.0°C	*** TEST PASSED ***
Control 2 IR	0.099%	12:46S	34.0°C	*** TEST PASSED ***
Ambient Air Blank	0.000%	12:47S		
Control 3 EC	0.099%	12:47S	34.0°C	*** TEST PASSED ***
Control 3 IR	0.098%	12:47S	34.0°C	*** TEST PASSED ***
Ambient Air Blank	0.000%	12:48S		

All tests within acceptable tolerance.

On this date, I installed the above indicated "NEW SOLUTION" in accordance with Alcotest 7110 operator training and procedures established by the (NJSP) Chief Forensic Scientist.

TEMPERATURE PROBE SERIAL NUMBER: DDWTP2-293_{gld}

Changed By:

Last Name: DELIANOCE

First Name: JOSEPH

MI: S

Signature: _____

Tr II J Delanoce # 6027

Badge No.: 6027

Date: 01/05/2011



State of New Jersey
 OFFICE OF THE ATTORNEY GENERAL
 DEPARTMENT OF LAW AND PUBLIC SAFETY
 DIVISION OF STATE POLICE
 POST OFFICE BOX 7068
 WEST TRENTON NJ 08628-0068
 (609) 882-2000

CHRIS CHRISTIE
 Governor

PAULA T. DOW
 Attorney General

KIM GUADAGNO
 Lt. Governor

COLONEL JOSEPH R. FUENTES
 Superintendent

CERTIFICATION OF ANALYSIS
0.10 PERCENT BREATH ALCOHOL SIMULATOR SOLUTION

ACCEPTANCE SPECIFICATIONS FOR BREATH ALCOHOL SIMULATOR SOLUTION: Ethyl alcohol concentration within, but not exceeding, the range of 0.117 to 0.125 grams per 100 milliliters of solution.

MANUFACTURER: Dräger Safety, Inc.


ANALYSIS DATE: 7/7/2010

BREATH ALCOHOL SIMULATOR SOLUTION LOT NUMBER: 10F079

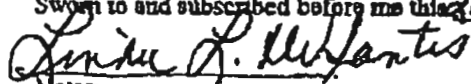
Representative samples of the above-referenced Lot Number were tested by Gas Chromatography and found to have a mean ethyl alcohol concentration range of 0.1182 to 0.1190 grams per 100 milliliters of solution.

This lot of breath alcohol simulator solution may be utilized as a known traceable standard for the purpose of conducting periodic tests, pursuant to N.J.A.C. 13:51-3.4, of approved breath test instruments (N.J.A.C. 13:51-3.5) utilized by law enforcement agencies in this State. The manufacturer's expiration date for this lot of breath alcohol simulator solution is June 10, 2012.

As Assistant Chief Forensic Scientist of the Division of State Police, I hereby certify and attest that the tests and results documented in this Certificate of Analysis were performed at my direction and under my supervision by personnel of, and at, the Office of Forensic Sciences of the Division of State Police on properly functioning and calibrated instruments and equipment. All procedures utilized are accurate, objective, and performed on a routine basis by personnel within the Office of Forensic Sciences, in accordance with their professional duties and responsibilities.


 Kenneth W. Kawalek, M.S.
 Assistant Chief Forensic Scientist
 Division of State Police

Sworn to and subscribed before me this 20th day of July, 2010.


 Notary

Linda L. Desantis
 Notary Public, New Jersey
 My Commission Expires 8-17-14



New Jersey Is An Equal Opportunity Employer
 Printed on Recycled Paper and Recyclable

